



ACCUJECT™

INJECTION SYSTEM

Your high-end solution



Medicel AG

Dornierstrasse 11

CH-9423 Altenrhein

SWITZERLAND


Tel. +41 71 727 10 50

Fax +41 71 727 10 55

info@medicel.com

www.medicel.com

medicel
SWISS TECHNOLOGY FOR SURGERY

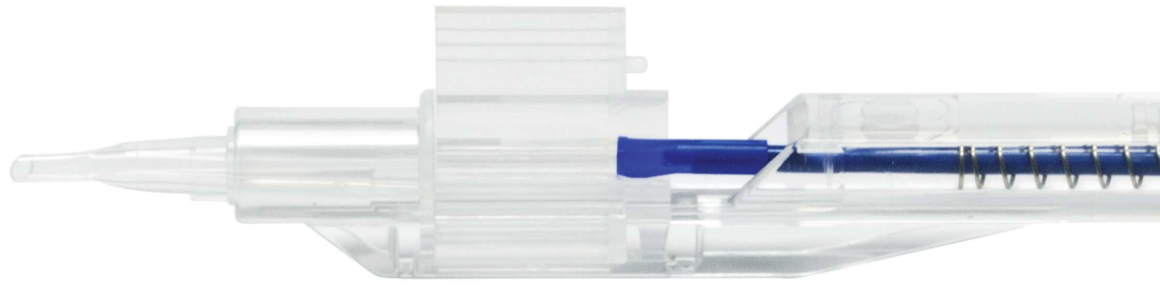

ACCUJECT™ / ACCUJECT™ SCREW


Medicel's ACCUJECT™ model is the most versatile injection system of all time. Its technologically sophisticated modular system enables maximum flexibility and adaptation to almost any lens design and every surgical technique. ACCUJECT™ allows the injection of 1- and 3-piece lenses, the use of hard- and soft-tips and loading the lens from the top (top-loaded) or from the back (back-loaded). Furthermore, it is also available as a syringe- or screw-type injector, and as a preloaded or nonpreloaded system for hydrophilic or hydrophobic lenses.

Incision sizes range from 3.2mm down to 1.6mm.

These countless versions are all based on a proven basic technology, which injects lenses repeatedly safe into the eye.


ACCUJECT™ - the all-rounder among lens injection systems.

Top-loaded version

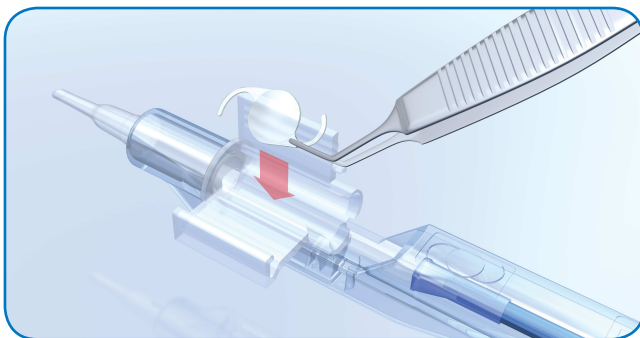

- Applicable for hydrophilic and hydrophobic lenses
- Large loading chamber
- Unrestricted access for IOL loading
- No force required for loading IOL
- Available as syringe- or screw-type
- Soft-tip and hard-tip plunger models available
- For 1-piece and 3-piece lenses
- Incision sizes 3.2mm to 1.6mm

Back-loaded version

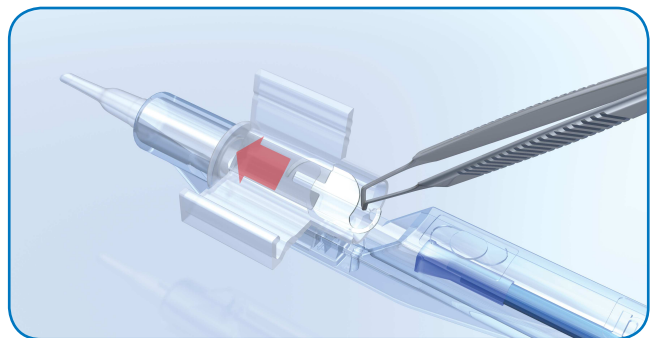

- Applicable for hydrophobic and some hydrophilic lenses
- Smart haptic folding for c-loop IOLs
- Available as syringe- or screw-type
- With proven soft-tip technique
- For 1-piece lenses
- Incision sizes 3.2mm to 1.8mm



ACCUJECT™ for 1-piece lenses



- Soft-tip solution (plunger) for 1-piece lenses

ACCUJECT™ for 3-piece lenses



- Hard-tip solution (plunger) for 3-piece lenses
- Haptic straightening mechanism

ACCUJECT™ PRO - preloaded system

ACCUJECT™ PRO for hydrophobic lenses



- For preloading hydrophobic (dry-stored) lenses
- Available as syringe- or screw-type
- Hard-tip version available upon request
- Smart haptic folding for c-loop IOLs
- Injector fully approved by Medicel
- Incision sizes from 3.2mm to 1.8mm

ACCUJECT™ PRO for hydrophilic lenses



- IOLs preloaded within the loading chamber
- Only loading chamber with IOL stored in liquid
- Available as syringe- or screw-type
- Low investment for IOL manufacturer
- Injector fully approved by Medicel
- Incision sizes from 3.2mm to 1.6mm

ACCUJECT™ - top-loaded

ITEM CODE	PRODUCT NAME	INCISION SIZE (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
LP604580	ACCUJECT™ 3.0-3P Injector Set	2.4	2.8	3.2	3-piece lenses
LP604520	ACCUJECT™ 2.5-3P Injector Set	2.2	2.5	2.8 - 2.9	3-piece lenses
LP604570	ACCUJECT™ 3.0-1P Injector Set	2.4	2.7-2.8	3.1	1-piece lenses
LP604500	ACCUJECT™ 2.6-1P Injector Set	2.2	2.5	2.8	1-piece lenses
LP604600	ACCUJECT™ 2.4-1P Injector Set	2.1	2.4	2.7	1-piece lenses
LP604590	ACCUJECT™ 2.1-1P Injector Set	1.9	2.1 - 2.2	2.4	1-piece lenses
LP604530	ACCUJECT™ 2.2-1P Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
LP604510	ACCUJECT™ 2.0-1P Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses
LP604540	ACCUJECT™ 1.8-1P Injector Set	1.8	2.0	2.2 - 2.3	1-piece lenses*
LP604550	ACCUJECT™ 1.6-1P Injector Set	1.6	1.8	2.0	1-piece lenses*

ACCUJECT™ - back-loaded

ITEM CODE	PRODUCT NAME	INCISION SIZE (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
LP604575	ACCUJECT™ 3.0-BL Injector Set	2.4	2.7 - 2.8	3.1	1-piece lenses
LP604505	ACCUJECT™ 2.6-BL Injector Set	2.2	2.5	2.8	1-piece lenses
LP604605	ACCUJECT™ 2.4-BL Injector Set	2.1	2.4	2.7	1-piece lenses
LP604535	ACCUJECT™ 2.2-BL Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
LP604595	ACCUJECT™ 2.1-BL Injector Set	1.9	2.1 - 2.2	2.4	1-piece lenses
LP604515	ACCUJECT™ 2.0-BL Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses
LP604545	ACCUJECT™ 1.8-BL Injector Set	1.8	2.0	2.2 - 2.3	1-piece lenses*

ACCUJECT™ SCREW - top-loaded

ITEM CODE	PRODUCT NAME	INCISION SIZE (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
AST3000	ACCUJECT™ SCREW 3.0-TL Injector Set	2.4	2.7 - 2.8	3.1	1-piece lenses
AST2400	ACCUJECT™ SCREW 2.4-TL Injector Set	2.1	2.4	2.7	1-piece lenses
AST2200	ACCUJECT™ SCREW 2.2-TL Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
AST2100	ACCUJECT™ SCREW 2.1-TL Injector Set	1.9	2.1 - 2.2	2.4	1-piece lenses
AST2000	ACCUJECT™ SCREW 2.0-TL Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses
AST1800	ACCUJECT™ SCREW 1.8-TL Injector Set	1.8	2.0	2.2 - 2.3	1-piece lenses*
AST1600	ACCUJECT™ SCREW 1.6-TL Injector Set	1.6	1.8	2.0	1-piece lenses*

ACCUJECT™ SCREW - back-loaded

ITEM CODE	PRODUCT NAME	INCISION SIZE (mm)			IOL TYPE
		wound-assisted	into-the-wound	into-the-bag	
ASB3000	ACCUJECT™ SCREW 3.0-BL Injector Set	2.4	2.7 - 2.8	3.1	1-piece lenses
ASB2600	ACCUJECT™ SCREW 2.6-BL Injector Set	2.2	2.5	2.8	1-piece lenses
ASB2400	ACCUJECT™ SCREW 2.4-BL Injector Set	2.1	2.4	2.7	1-piece lenses
ASB2200	ACCUJECT™ SCREW 2.2-BL Injector Set	1.9 - 2.0	2.2	2.5	1-piece lenses
ASB2100	ACCUJECT™ SCREW 2.1-BL Injector Set	1.9	2.1 - 2.2	2.4	1-piece lenses
ASB2000	ACCUJECT™ SCREW 2.0-BL Injector Set	1.8 - 1.9	2.1	2.3 - 2.4	1-piece lenses

*hydrophilic only