



POSTERIOR SURGICAL DYES: FOCUS ON PURITY

MEMBRANEBLUE-DUAL® and ILM-BLUE®

Pure performance: Focus on ILM-BLUE® and MEMBRANEBLUE-DUAL®

- Next generation posterior dyes.
- Optimal injection characteristics.
- ILM-BLUE® for selective ILM staining.
- MEMBRANEBLUE-DUAL® for effective ILM and ERM staining.
- High purification for patient safety.
- Ready to use in syringes.

ILM-BLUE®

- ILM-BLUE® is a dye to stain the ILM in vitreoretinal surgery.
- By injecting ILM-BLUE® in the vitreous cavity the ILM will be clearly stained and easily distinguished from the underlying, unstained retina and can be removed selectively.
- Due to an integrated carrier 4% PEG solution, ILM-BLUE® can be injected in a BSS filled eye and sinks immediately as a cohesive ball to the fundus of the eye and only stains the targeted tissue without diffusion throughout the whole globe.

Specifications:


Composition of one 0.5ml syringe:

- 0.125 mg Brilliant Blue G
- PEG 4% PEG 3350

Concentration: 0.25 g/l
pH-value: 7.3 – 7.6
Osmolality: 301 - 369 mOsm/kg H₂O

ILM-BLUE®: Posterior Dye for staining the ILM consists of:

- BBG 0,025% + 4% PEG
- Highly purified BBG up to 97% purity!

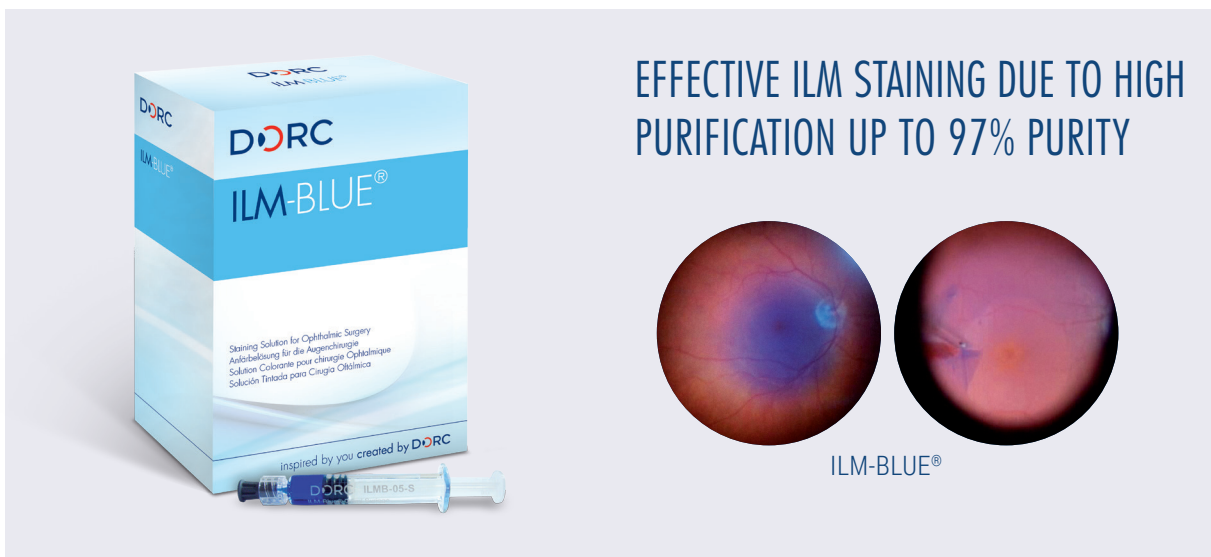

ILM.05S

ILM-BLUE®, 0.5ml Syringe (Sterile, box/5).

SAFETY TESTS

- Toxicity tests on ARPE-19 human retinal pigment epithelium cells showed an advantageous safety profile of MEMBRANEBLUE-DUAL® and ILM-BLUE®.
- Increased purity of the dyes as well as the addition of PEG3350 results in even safer products.
- Safety testing was complemented by electron microscopy ultrastructural analysis and in vitro cytotoxicity tests according to ISO-10993.

Advanced injection CHARACTERISTICS due to "carrier" 4% PEG SOLUTION

- Creates better injection characteristics under BSS.
- Sinks immediately to posterior pole as a cohesive ball.
- Stains only target tissue and no diffusion throughout the eye.

MEMBRANEBLUE-DUAL®

- MEMBRANEBLUE-DUAL® is a dye to stain both the ILM as well as the ERM and PVR membrane, without compromising the staining effect within one injection.
- By injecting MEMBRANEBLUE-DUAL® in the vitreous cavity the membrane will be clearly stained and easily distinguished from the underlying, unstained retina and can be removed selectively.
- MEMBRANEBLUE-DUAL® is a stable mix and stains the ILM at the same level as ILM-Blue® and stains the ERM and PVR at a higher level than the classic MembraneBlue.
- Due to an integrated carrier 4% PEG solution, MEMBRANEBLUE-DUAL® can be injected in a BSS filled eye and sinks immediately as a cohesive ball to the fundus of the eye and only stains the targeted tissue without diffusion throughout the whole globe.

Specifications:

Composition of one 0.5ml syringe:

- 0.125 mg Brilliant Blue G
- 0.75 mg Trypan Blue
- PEG 4% PEG 3350

Concentration: 1.75 g/l
pH-value: 7.3 – 7.6
Osmolality: 301-369 mOsm/kg H2O

MEMBRANEBLUE-DUAL®: Dye for staining ILM, ERM & PVR membranes consists of:

- Combination of TrypanBlue 0,15% + BBG 0,025% + 4% PEG
- Highly purified TB + BBG up to 97%



MBD.05S

MEMBRANEBLUE-DUAL® 0.5 ml Syringe (Sterile, box/5).

PURITY & SAFETY

- ILM-BLUE® & MEMBRANEBLUE-DUAL® are produced according to Current Good Manufacturing Practice (CGMP)
- ILM-BLUE® & MEMBRANEBLUE-DUAL® are produced with highly purified BBG + Trypanblue at a drug substance level.

PATENT PROTECTED

USA 6,696,430B1
6,372,449
20080090914A1
Australia 75731 8
Japan 4200222
Europe EP 1 819 366 B1



Technology
Licensed from
AQUMEN, Japan

In Austria please contact:
D.O.R.C. GmbH
Phone: +43 1 59 99 96 30
E-mail: info-at@dorc.eu

In Brasil please contact:
D.O.R.C. do Brasil Produtos e
Serviços Oftalmológicos Ltda
Phone: +55 11 22 93 36 05
E-Mail: vendas@dorc.com.br

In France please contact:
D.O.R.C. France S.A.R.L.
Phone: +33 1 70 82 47 67
E-mail: ventes@dorc.eu

In Germany please contact:
D.O.R.C. Deutschland GmbH
Phone: +49 211 52 29 15 80
E-mail: bestellung@dorc.eu


In Italy please contact:
D. O. R. C. Italy S. r. l.
Phone: +39 02 00 69 64 01
E-mail: ordini-it@dorc.eu

In Scandinavia please contact:
D.O.R.C. Scandinavia AB
Phone: +46 19 12 15 14
E-mail: order@dorc.eu

In Spain please contact:
D.O.R.C. España, S. L.
Phone: +34 900 64 90 06 or
+34 915 03 00 13
E-Mail: pedidos-es@dorc.eu

In UK please contact:
D.O.R.C. Limited
Phone: +44 80 00 51 83 77
E-mail: orders-uk@dorc.eu

AUTHORIZED D.O.R.C. DEALER



Scan & download
the D.O.R.C. app



For additional information please contact:
D.O.R.C. Dutch Ophthalmic Research Center (International) B.V.
Scheijdelveweg 2, 3214 VN Zuidland, The Netherlands

Phone: +31 181 45 80 80
Fax: +31 181 45 80 90
E-mail: sales@dorc.eu